


Rolyan®

CuraDrape™

Performance Health are pleased to announce the introduction of a new thermoplastic splinting material to join the Rolyan range.

CuraDrape is easy to cut, easy to form and lightweight with a velvety soft, silky smooth finish. With a minimum to moderate resistance to stretch **CuraDrape** is highly conformable, has excellent drapability and requires minimal handling.

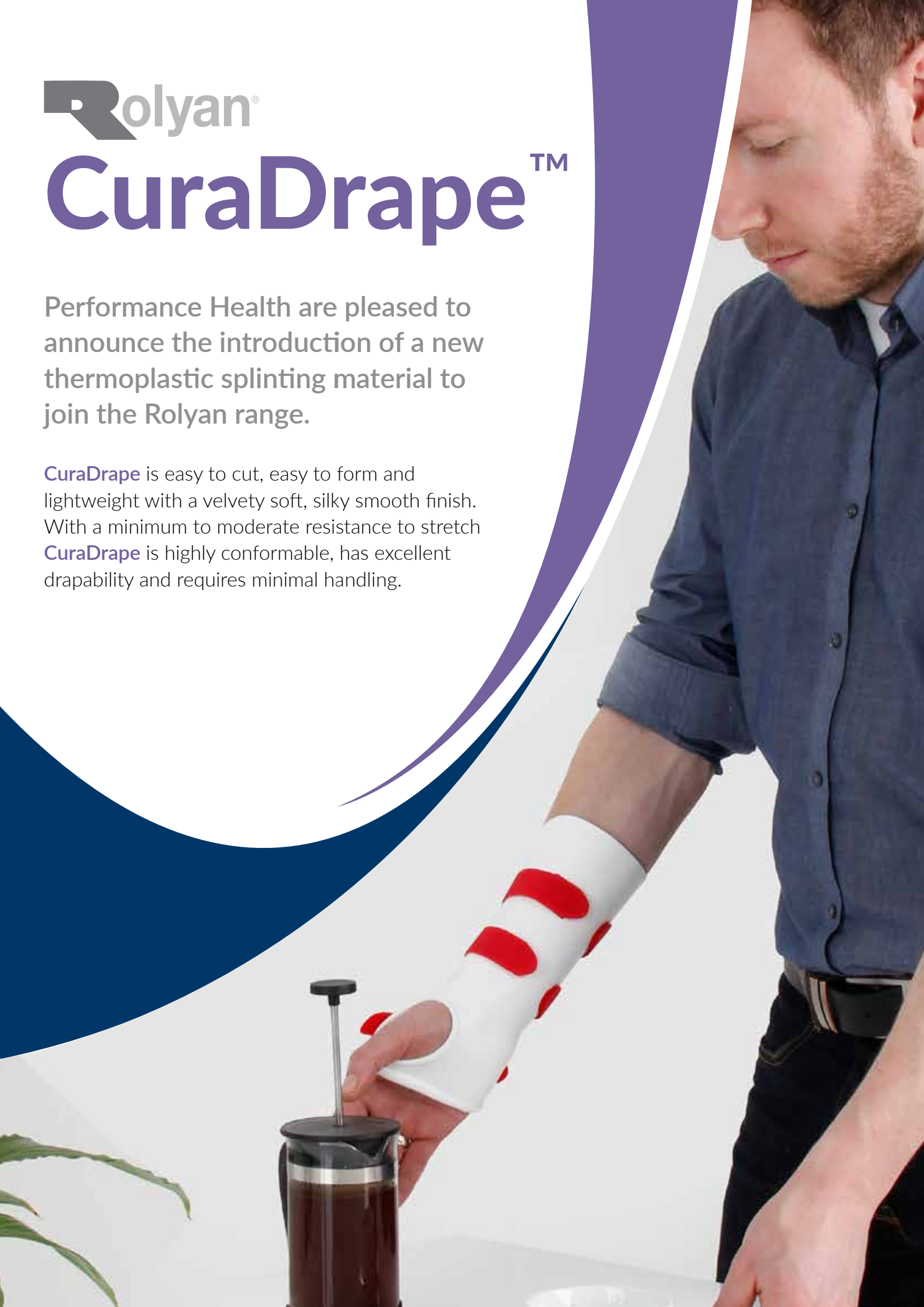



CuraDrape™

Features

- Lightweight
- Highly drapeable
- Minimum to moderate resistance to stretch
- Brilliant bright white and smooth finish
- Moderate memory
- Temporary bond



Thickness
Available in 1.6mm, 2.4mm, 3.2mm

Perforations
Solid, 1% perf, Optiperf and Ultraparf

Colour
White

Choice


With the broadest range of high-quality splinting materials in the industry, Performance Health makes it easy to find a product for every splint, every therapist and every experience level.

Support

From our large, knowledgeable field sales representatives through to customer services and clinical support, with Performance Health you'll find the expertise, and information you need to choose the right product for every patient, every time.

Quality

Rolyan® Manufacturing Company began as a small family business in 1967. All our thermoplastic materials continue to be manufactured in the USA to a high quality and ethical standards. Rolyan®'s extensive range of thermoplastics has continued to grow and develop over the years.

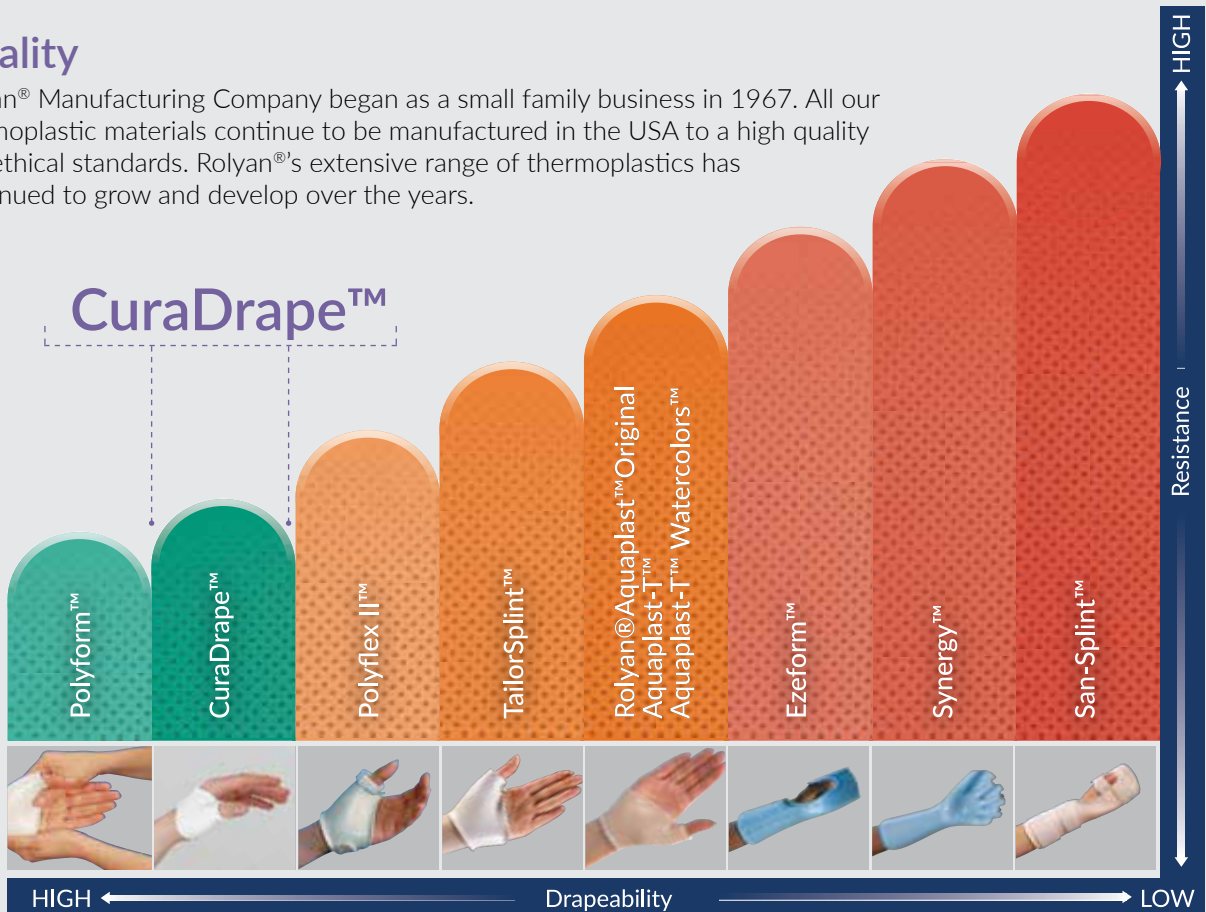



Rolyan® CuraDrape™

Minimum Resistance to Stretch / Coated / Latex Free / Made in USA



1.6mm 65°C to 70°C for 30 seconds
2.4mm 65°C to 70°C for 30 seconds
3.2mm 65°C to 70°C for 1 minute



1.6mm 1 minute of working time
2.4mm 1 to 2 minutes of working time
3.2mm 3 to 5 minutes working time



Handling Characteristics

- Excellent drapability, conforming well to contours providing an intimate fit
- Requires minimal handling
- Minimum to moderate stretch
- Minimal handling required
- Easy trim edges that finish well
- Coated: allowing temporary bonds

Bond-Ability

- Coated: allowing temporary bonds for a temporary bond pinch heated surface together
- For a permanent bond remove the coating with spirits or by safely scarping the surface with a sharp instrument

Clinical Recommendations

Ideal for conditions that cannot tolerate firm handling or require a lightweight material ie:


- Acute trauma
- Pain and joint pathologies
- Thumb, finger, wrist and hand based splints
- Paediatrics

Ordering information

Thickness	Sheet Qty	Colour	Approx. Size	Code By Perforations			
				Solid	1% Perf	13% UltraPerf	19% OptiPerf
1.6mm	1	White	30 x 46cm			PAT-081707793	
2.4mm	1	White	46 x 61cm	PAT-081707538	PAT-081707553		
3.2mm	1	White	46 x 61cm	PAT-081707579	PAT-081707595		PAT-081707611



For samples of CuraDrape™
or any Rolyan Thermoplastics,
contact your Performance Health
Area Sales Manager now!

 1300 473 422

 sales@performhealth.com.au

 www.performhealth.com.au



PERFORMANCE
HEALTH®

Partner of:



australianhand
therapyassociation^{inc.}